



DICOM Conformance Statement

I-SoftView 6.3

Document Number: 2004-02987

Revision: 5.0

March 24, 2009

I-SoftView is sold by Toshiba Medical Systems Europe
under the marketing name ReportDirect



DICOM Conformance Statement for Cedara I-SoftView™ 6.3

Document Number: 2004-02987

Revision: 5.0

March 24, 2009

Copyright © 2009 by Merge Healthcare

(This page has been left blank intentionally.)

Table of Contents

1. INTRODUCTION.....5

1.1 PURPOSE OF THIS DOCUMENT.....5

1.2 RELATED DOCUMENTS5

1.3 DEFINITIONS5

1.4 ACRONYMS AND ABBREVIATIONS5

1.5 IMPORTANT NOTE TO THE READER5

2. IMPLEMENTATION MODEL.....7

2.1 APPLICATION DATA FLOW DIAGRAM7

2.2 FUNCTIONAL DEFINITIONS OF APPLICATION ENTITIES9

2.3 SEQUENCING OF REAL-WORLD ACTIVITIES10

3. APPLICATION ENTITY SPECIFICATIONS.....11

3.1 DCSERVER APPLICATION ENTITY SPECIFICATION11

3.1.1 *Association Establishment Policies*12

3.1.1.1 General.....12

3.1.1.2 Number of Associations.....12

3.1.1.3 Asynchronous Nature.....13

3.1.1.4 Implementation Identifying Information.....13

3.1.2 *Association Initiation Policy*.....13

3.1.2.1 Query Request.....13

3.1.2.1.1 Associated Real World Activity13

3.1.2.1.2 Proposed Presentation Contexts13

3.1.2.2 Move Request16

3.1.2.2.1 Associated Real World Activity16

3.1.2.2.2 Proposed Presentation Contexts16

3.1.2.3 Store Request.....17

3.1.2.3.1 Associated Real World Activity17

3.1.2.3.2 Proposed Presentation Contexts17

3.1.3 *Association Acceptance Policy*.....20

3.1.3.1 Storage Association Request.....20

3.1.3.1.1 Associated Real-World Activity20

3.1.3.1.2 Presentation Context Table.....20

3.1.3.1.3 Presentation Context Acceptance Criterion27

3.1.3.1.4 Transfer Syntax Selection Policies27

3.1.3.2 Query Association Request.....27

3.1.3.2.1 Associated Real-World Activity27

3.1.3.2.2 Presentation Context Table.....27

3.1.3.2.3 Presentation Context Acceptance Criterion29

3.1.3.2.4 Transfer Syntax Selection Policies29

3.1.3.3 Move Association Request.....29

3.1.3.3.1 Associated Real-World Activity29

3.1.3.3.2 Presentation Context Table.....29

3.1.3.3.3 Presentation Context Acceptance Criterion30

3.1.3.3.4 Transfer Syntax Selection Policies30

3.2 HCSERVER APPLICATION ENTITY SPECIFICATION31

3.2.1 *Association Establishment Policies*31

3.2.1.1 General.....31

3.2.1.2 Number of Associations.....31

3.2.1.3 Asynchronous Nature.....31

3.2.1.4 Implementation Identifying Information.....31

3.2.2 *Association Initiation Policy*.....31

3.2.2.1 Print to remote printer.....31

3.2.2.1.1 Associated Real World Activity31

3.2.2.1.2 Proposed Presentation Contexts 31

3.2.3 Association Acceptance Policy 35

4. COMMUNICATION PROFILES 36

4.1 SUPPORTED COMMUNICATION STACKS (PARTS 8,9)..... 36

4.1.1 TCP/IP Stack 36

4.1.1.1 API..... 36

4.1.1.2 Physical Media Support 36

5. EXTENSIONS/SPECIALIZATION'S/PRIVATIZATION'S 37

5.1 STANDARD/EXTENDED/SPECIALIZED/PRIVATE SOPS 37

5.2 PRIVATE TRANSFER SYNTAX'S 37

6. CONFIGURATION 38

6.1 DCSERVER APPLICATION ENTITY CONFIGURATION 38

6.2 HCSERVER APPLICATION ENTITY CONFIGURATION 38

7. SUPPORT OF EXTENDED CHARACTER SETS 39

1. Introduction

1.1 Purpose of this Document

This document is the DICOM Conformance Statement for Cedara I-SoftView™ 6.3. The purpose of this document is to describe how the workstation interacts with other DICOM devices on the network.

1.2 Related Documents

The Digital Imaging and Communications in Medicine (DICOM) standard. NEMA PS 3.1-13 2003 and Supplements.

1.3 Definitions

DICOM terms are used throughout this Conformance Statement. For a description of these, consult the DICOM standard publication.

Word	Definition
dcserver	The executable name of the Cedara DICOM Image Transfer Server
hcserver	The executable name of the Cedara Hardcopy Server

1.4 Acronyms and Abbreviations

DICOM abbreviations are used throughout this Conformance Statement. For a description of these, consult the DICOM standard publication.

Acronym	Meaning

1.5 Important note to the reader

The use of this conformance statement by itself does not guarantee successful interoperability of Merge OEM products with equipment from other vendors. The user or integrator of Merge OEM products should keep the following issues in mind:

1. Successful interoperability of the workstation with other devices may require functions that are not specified within the scope of DICOM. It is the user's or integrator's responsibility to ensure that the proper analysis and validation is performed to verify the connection.
2. Test procedures should be used to verify that data transferred into the workstation is correct. This is often done with the aid of phantom scans using a variety of clinical protocols.
3. Test procedures should be used to verify connectivity. Issues such as full database and broken connections should be verified.
4. The current version of this DICOM Conformance Statement is also available on the Merge OEM web page at <http://www.MergeOEM.com/>

5. The DICOM standard will continually evolve to meet new user requirements. Merge OEM reserves the right to make changes to its products or to discontinue its delivery. The user or integrator should ensure that any non-Merge OEM device providers, which connect with Merge OEM devices, should also follow the standard. Failure to do so will likely result in future connectivity problems.

2. Implementation Model

The Cedara I-SoftView 6.3 from Merge OEM is a medical diagnostic and review station. The workstation uses DICOM services to:


- import and export images
- query the content of other devices and initiate transfers
- print images

The DICOM services are implemented using two Application Entities. The DICOM SERVER AE is used for image transfers and is implemented by a Windows application called *dcserver*. The HARDCOPY SERVER AE is used for printing and is implemented by a Windows application called *hcserver*.

2.1 Application Data Flow Diagram

The diagram for the Implementation Model is shown in Figure 1.

Figure 1. Application Entity Implementation Model

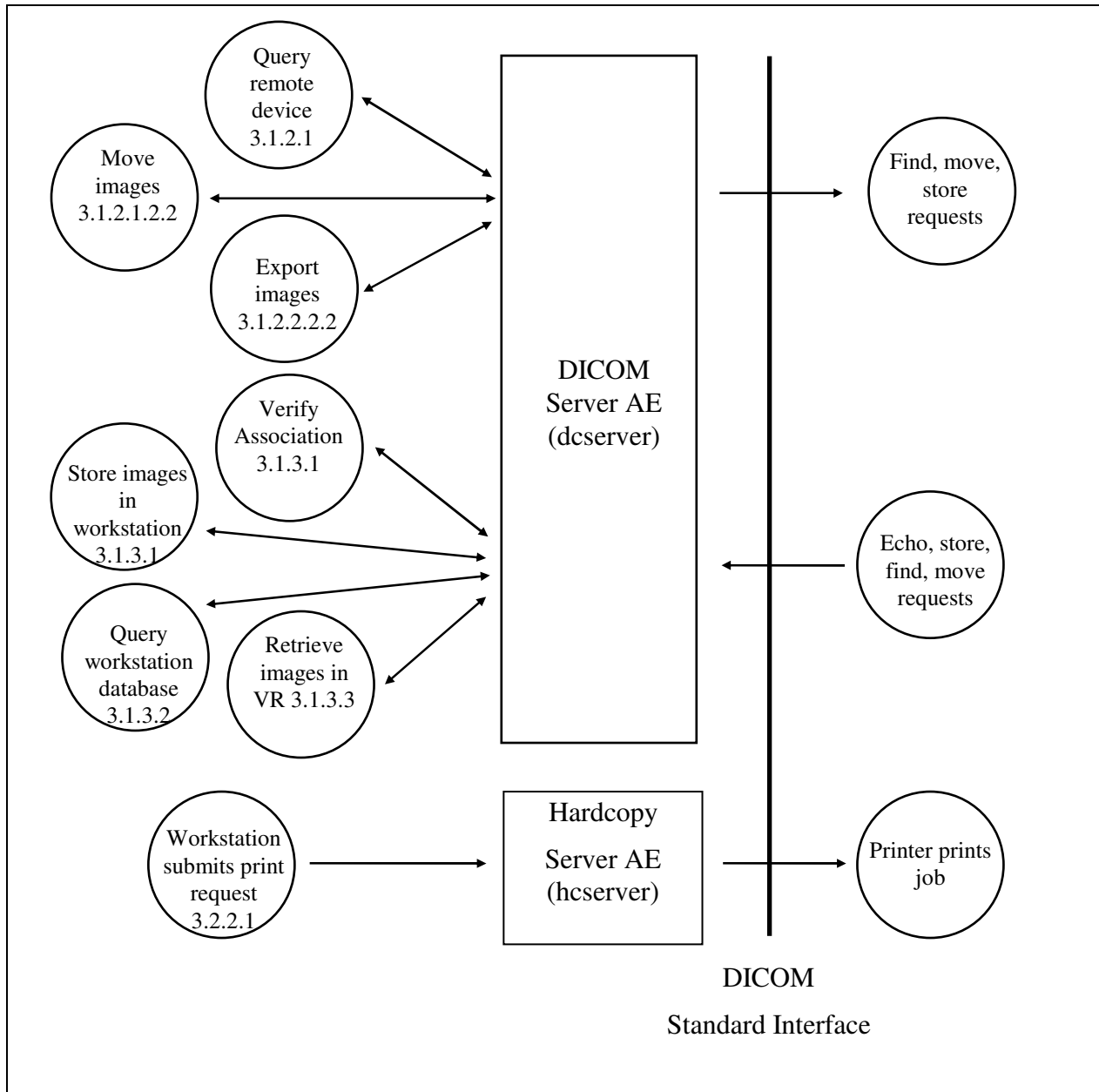

The *dcserver* is expected to be running on the local workstation. A Remote Application Entity initiates an association for Storage Services. Upon notification of acceptance of the association parameters, the Remote Application Entity sends Information Objects to the *dcserver* that stores them in a local database for future use by the workstation.

The user initiates Query and Retrieve requests using the *dcserver* component, generally by interaction with a User Interface. The *dcserver* component initiates an association with the Remote Application Entity and uses the Query or Retrieve Service Classes to issue commands. The Remote Application Entity responds as a Query/Retrieve Service Class Provider performing C-FIND and C-MOVE operations as required. The *dcserver* component passes the status responses for these commands to the workstation User Interface for interpretation and display.

The *hcserver* can be running on the local or a remote workstation. The application using the *hcserver* requests printing to a print device. The *hcserver* initiates an association with a DICOM print SCP for the purpose of printing the job requested by the application. The *hcserver* can handle simultaneous associations with a number of DICOM print SCPs.

2.2 Functional Definitions of Application Entities

The *dcserver* component operates as a daemon. The start-up sequence of the system initiates its execution. The *dcserver* is left running whether the workstation is operational or not.

The *dcserver* uses a configuration file that contains information used to validate association attempts from Remote Application entities. The *dcserver* then listens on the configured port for association requests.

An association request for Storage Services from a Remote Application Entity causes *dcserver* to validate the request according to the configuration parameters set at execution time. The Remote Application Entity then sends the Information Object Instance. The *dcserver* stores the received Information Object Instance in its local database if the data does not already exist. The data remains in the database until removed by some action external to this Application Entity.

An association request from a Remote Application Entity for Query or Move Services causes *dcserver* to validate the request according to the configuration parameters set at execution time. The Remote Application Entity then sends the Query or Retrieve request. The *dcserver* searches the local database for the instance(s) specified. If the request was C_FIND, then a response is returned for each match. If the request was C_MOVE, then an association is originated to the destination Application Entity specified in the C_MOVE message. Incremental responses are sent to the C_MOVE originator to indicate progress of the request.

A request from the workstation User Interface causes the *dcserver* component to initiate an association with a Remote Application Entity. The Service Classes offered are specified in the configuration file. The user can then initiate query and retrieve requests to *dcserver* that are sent to the Remote Application Entity. The workstation User Interface displays the responses from the Remote Application Entity.

The *hcserver* component operates as a daemon. The start-up sequence of the system initiates its execution. The *hcserver* is left running whether the workstation is operational or not.

The *hcserver* uses a configuration file to determine the list of printer devices connected to the server and the properties of each printer.

Association and release requests are logged to the Windows event log or a log file on Windows 95. Various error and warning indications are also logged using *syslogd*.

2.3 Sequencing of Real-World Activities

Not applicable.

3. Application Entity Specifications

3.1 Dcserver Application Entity Specification

The workstation DICOM Image Transfer capability consists of two logical components. The SCU portion originates associations for Store, Query and Retrieve operations. The SCP portion accepts associations for Store, Query and Retrieve operations. The two components are configured with the same Application Entity Title for use in the workstation. They are treated as a single Application Entity in this description.

The *dcserver* Application Entity provides Standard Conformance to the following DICOM V3.1 SOP Classes as an SCP:

SOP Class Name	SOP Class UID
Verification	1.2.840.10008.1.1
Patient Root Query/Retrieve Model - FIND	1.2.840.10008.5.1.4.1.2.1.1
Patient Root Query/Retrieve Model - MOVE	1.2.840.10008.5.1.4.1.2.1.2
Study Root Query/Retrieve Model - FIND	1.2.840.10008.5.1.4.1.2.2.1
Study Root Query/Retrieve Model - MOVE	1.2.840.10008.5.1.4.1.2.2.2
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1
Computed Tomography Image Storage	1.2.840.10008.5.1.4.1.1.2
Digital Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1.1
Mammography Image Storage	1.2.840.10008.5.1.4.1.1.1.2
Magnetic Resonance Image Storage	1.2.840.10008.5.1.4.1.1.4
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20
Nuclear Medicine Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.5
Positron Emission Tomography Image Storage	1.2.840.10008.5.1.4.1.1.128
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7
Ultrasound Image Storage	1.2.840.10008.5.1.4.1.1.6.1
Ultrasound Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.6
Ultrasound Multi-frame Image Storage	1.2.840.10008.5.1.4.1.1.3.1
Ultrasound Multi-frame Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.3
X-Ray Angiographic Bi-Plane Image Storage	1.2.840.10008.5.1.4.1.1.12.3
X-Ray Angiographic Image Storage	1.2.840.10008.5.1.4.1.1.12.1
X-Ray Radiofluoroscopic Image Storage	1.2.840.10008.5.1.4.1.1.12.2
Greyscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1
Key Image Note	1.2.840.10008.5.1.4.1.1.88.59
Mammography CAD SR IOD	1.2.840.10008.5.1.4.1.1.88.50

The *dcs* Application Entity provides Standard Conformance to the following DICOM V3.1 SOP Classes as an SCU:

SOP Class Name	SOP Class UID
Patient Root Query/Retrieve Model - FIND	1.2.840.10008.5.1.4.1.2.1.1
Patient Root Query/Retrieve Model - MOVE	1.2.840.10008.5.1.4.1.2.1.2
Study Root Query/Retrieve Model - FIND	1.2.840.10008.5.1.4.1.2.2.1
Study Root Query/Retrieve Model - MOVE	1.2.840.10008.5.1.4.1.2.2.2
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1
Computed Tomography Image Storage	1.2.840.10008.5.1.4.1.1.2
Digital Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1.1
Mammography Image Storage	1.2.840.10008.5.1.4.1.1.1.1
Magnetic Resonance Image Storage	1.2.840.10008.5.1.4.1.1.4
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20
Nuclear Medicine Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.5
Positron Emission Tomography Image Storage	1.2.840.10008.5.1.4.1.1.128
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7
Ultrasound Image Storage	1.2.840.10008.5.1.4.1.1.6.1
Ultrasound Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.6
Ultrasound Multi-frame Image Storage	1.2.840.10008.5.1.4.1.1.3.1
Ultrasound Multi-frame Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.3
X-Ray Angiographic Bi-Plane Image Storage	1.2.840.10008.5.1.4.1.1.12.3
X-Ray Angiographic Image Storage	1.2.840.10008.5.1.4.1.1.12.1
X-Ray Radiofluoroscopic Image Storage	1.2.840.10008.5.1.4.1.1.12.2
Greyscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1
Key Image Note	1.2.840.10008.5.1.4.1.1.88.59
Mammography CAD SR IOD	1.2.840.10008.5.1.4.1.1.88.50

3.1.1 Association Establishment Policies

3.1.1.1 General

The User of the workstation User Interface can select which Application Entity to associate with for Query and Retrieve operations. The configuration file contains the configuration parameters such as host name, port number and specific SOP Classes to negotiate for each accessible Application Entity.

The maximum PDU size is 65536 bytes.

The *dcs* Application Entity always proposes or accepts the Verification SOP Class.

3.1.1.2 Number of Associations

The *dcs* can initiate multiple associations concurrently.

A configuration parameter is provided to limit the number of associations that can be originated. The default is 3.

Another configuration parameter is provided to limit the number of associations that *dcserver* can accept. The default is 10.

3.1.1.3 Asynchronous Nature

The application does not support asynchronous operations and will not perform asynchronous window negotiation.

3.1.1.4 Implementation Identifying Information

The *dcserver* implementation class UID is 2.16.124.113531.1.1.

The *dcserver* implementation version name is CEDARA VR 4.0.

3.1.2 Association Initiation Policy

This section details the action of the *dcserver* SCU component as a result of user initiated activity on the workstation User Interface.

3.1.2.1 Query Request

3.1.2.1.1 Associated Real World Activity

The user of the workstation selects the Query operation button on the user interface. The user can specify wild card or specific information for Patient Name, Patient ID, Study ID, Study Date range or Referring Physician.

Wild card queries can result in excessive number of responses. The user interface is able to restrict the number of patients displayed. A warning dialogue is displayed to indicate that too many matches were found. A query cancel command is sent when the limit is reached.

The DICOM transfer utility defaults to using Study Root Query Model when initiating query request. The query model used can be changed to Patient Root Query Model by changing a configuration parameter.

3.1.2.1.2 Proposed Presentation Contexts

The following table describes the Presentation Contexts that may be presented for the Query request. The configuration file contains 1 of the listed Abstract Syntax's.

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Patient Root Query / Retrieve Model – FIND	1.2.840.10008.5.1.4.1.2.1.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Study Root Query / Retrieve Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

3.1.2.1.2.1 SOP Specific Conformance for Patient Root Query/Retrieve Model – FIND

The *dcserver* does not use Relational Queries.

The *dcserver* does not use Extended Negotiation.

The Keys supported are listed below:

Patient Level Keys

Description	Tag	Type
Patient’s Name*	(0010,0010)	R
Patient ID*	(0010,0020)	U
Patient’s Birth Date	(0010,0030)	O
Patient’s Birth Time	(0010,0032)	O
Patient’s Sex	(0010,0040)	O
Other Patient Ids	(0010,1000)	O
Other Patient Names	(0010,1001)	O
Ethnic Group	(0010,2160)	O
Patient Comments	(0010,4000)	O

Study Level Keys

Description	Tag	Type
Study Date*	(0008,0020)	R
Study Time*	(0008,0030)	R
Accession Number	(0008,0050)	R
Study ID*	(0020,0010)	R
Study Instance UID	(0020,000D)	U
Referring Physician’s Name*	(0008,0090)	O
Study Description	(0008,1030)	O
Name of Physician(s) Reading Study	(0008,1060)	O
Admitting Diagnoses Description	(0008,1080)	O
Patient’s Age	(0010,1010)	O
Patient’s Size	(0010,1020)	O
Patient’s Weight	(0010,1030)	O
Occupation	(0010,2180)	O
Additional Patient History	(0010,21B0)	O

Description	Tag	Type
Other Study Numbers	(0020,1070)	O
Interpretation Author	(4008,010C)	O

Series Level Keys

Description	Tag	Type
Modality*	(0008,0060)	R
Series Number*	(0020,0011)	R
Series Instance UID	(0020,000E)	U

Image Level Keys

Description	Tag	Type
Image Number	(0020,0013)	R
SOP Instance UID	(0008,0018)	U

- The keys marked with an asterisk are displayed in the DICOM transfer utility.

3.1.2.1.2.2 SOP Specific Conformance for Study Root Query/Retrieve Model – FIND

The *dcserver* does not use Relational Queries.

The *dcserver* does not use Extended Negotiation.

The Keys supported are listed below:

Study Level Keys

Description	Tag	Type
Study Date*	(0008,0020)	R
Study Time*	(0008,0030)	R
Accession Number	(0008,0050)	R
Patient's Name*	(0010,0010)	R
Patient ID*	(0010,0020)	R
Study ID*	(0020,0010)	R
Study Instance UID	(0020,000D)	U
Referring Physician's Name	(0008,0090)	O
Study Description	(0008,1030)	O
Name of Physician(s) Reading Study	(0008,1060)	O
Admitting Diagnoses Description	(0008,1080)	O
Patient's Birth Date	(0010,0030)	O
Patient's Birth Time	(0010,0032)	O
Patient's Sex	(0010,0040)	O
Other Patient Ids	(0010,1000)	O
Other Patient Names	(0010,1001)	O

Description	Tag	Type
Patient's Age	(0010,1010)	O
Patient's Size	(0010,1020)	O
Patient's Weight	(0010,1030)	O
Ethnic Group	(0010,2160)	O
Occupation	(0010,2180)	O
Additional Patient History	(0010,21B0)	O
Patient Comments	(0010,4000)	O
Other Study Numbers	(0020,1070)	O
Number of Study Related Images	(0020,1208)	O
Interpretation Author	(4008,010C)	O

Series Level Keys

Description	Tag	Type
Modality*	(0008,0060)	R
Series Number*	(0020,0011)	R
Series Instance UID	(0020,000E)	U

Image Level Keys

Description	Tag	Type
Image Number	(0020,0013)	R
SOP Instance UID	(0008,0018)	U

The keys marked with an asterisk are displayed in the DICOM transfer utility.

3.1.2.2 Move Request

3.1.2.2.1 Associated Real World Activity

The user selects one or more studies and/or series within studies from a list presented as a result of a previous Query operation.

The user of the workstation then selects the Send operation button on the user interface to initiate the move operation. The Destination Application Entity Title is selectable on the User Interface.

3.1.2.2.2 Proposed Presentation Contexts

The following table describes the Presentation Contexts that may be presented for the Move request. The configuration file contains 1 of the listed Abstract Syntax's.

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Patient Root Query / Retrieve Model – MOVE	1.2.840.10008.5.1.4.1.2.1.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Study Root Query / Retrieve Model – MOVE	1.2.840.10008.5.1.4.1.2.2.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

3.1.2.2.2.1 SOP Specific Conformance for Patient Root Query/Retrieve Model – MOVE

This implementation supports transfers against the Patient Query/Retrieve Information Model described in Section C.6.1.1 of NEMA PS3.4 (1996) Annex C using the C-MOVE SCU behaviour described in Section C.4.2.2 of NEMA PS3.4 (1996) Annex C.

3.1.2.2.2.2 SOP Specific Conformance for Study Root Query/Retrieve Model – MOVE

This implementation supports transfers against the Study Query/Retrieve Information Model described in Section C.6.2.1 of NEMA PS3.4 (1994) Annex C using the C-MOVE SCU behaviour described in Section C.4.2.2 of NEMA PS3.4 (1994) Annex C.

3.1.2.3 Store Request

3.1.2.3.1 Associated Real World Activity

The *dcserver* Application Entity initiates an association for C_STORE if it has received a valid C_MOVE message from a Remote Application Entity. The SOP Class UID of the Information Object to be sent over the C_STORE context is used to verify that a valid Presentation Context exists prior to issuing the C_STORE message. A mismatch results in no message being sent but the association remains active.

3.1.2.3.2 Proposed Presentation Contexts

The following table describes the Presentation Contexts that may be presented for the Store request. The configuration file contains 1 or more of the listed Abstract Syntax's.

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Computed Tomography Image Storage	1.2.840.10008.5.1.4.1.1.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Digital Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Mammography Image Storage	1.2.840.10008.5.1.4.1.1.1.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Magnetic Resonance Image Storage	1.2.840.10008.5.1.4.1.1.4	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Nuclear Medicine Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.5	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Positron Emission Tomography Image Storage	1.2.840.10008.5.1.4.1.1.128	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Ultrasound Image Storage	1.2.840.10008.5.1.4.1.1.6.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Ultrasound Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.6	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Ultrasound Multi-frame Image Storage	1.2.840.10008.5.1.4.1.1.3.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Ultrasound Multi-frame Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.3	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
X-Ray Angiographic Bi-Plane Image Storage	1.2.840.10008.5.1.4.1.1.12.3	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
X-Ray Angiographic Image Storage	1.2.840.10008.5.1.4.1.1.12.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
X-Ray Radiofluoroscopic Image Storage	1.2.840.10008.5.1.4.1.1.12.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Greyscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Key Image Note	1.2.840.10008.5.1.4.1.1.88.59	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None
Mammography CAD SR IOD	1.2.840.10008.5.1.4.1.1.88.50	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

3.1.2.3.2.1 SOP Specific Conformance for Image Storage SOP Classes

This implementation supports transfers as an SCU as described in NEMA PS3.4 Annex B.

The status returned by the accepting Application Entity is used to indicate success or failures of the C_MOVE sub-operation that initiated the transfer. In no case is the Information Object deleted from the local database.

Extended negotiation is not used by *dcserver* for this SOP Class.

The *dcserver* can be configured to export private Merge OEM data attributes. Group 29 is used for such attributes. The Private Data Creator element value is "Silhouette VRS 3.0"

3.1.3 Association Acceptance Policy

Parameters in the *dcserver* configuration file determine association acceptance. Association acceptance can be controlled on the basis of Called Application Entity Title, Calling Application Entity Title and SOP Class UID matching. Acceptance control ranges from no limitations to very specific acceptance policies.

A configuration parameter can be set to limit the number of accepted associations to a specific value.

3.1.3.1 Storage Association Request

3.1.3.1.1 Associated Real-World Activity

The *dcserver* stores image Information Object Instances received on the accepted association into its attached database.

3.1.3.1.2 Presentation Context Table

The following table lists the possible Presentation Contexts. The Application Entity configuration file specifies which of these Presentation Contexts are actually used in a specific configuration.

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Verification	1.2.840.10008.1.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Computed Tomography Image Storage	1.2.840.10008.5.1.4.1.1.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Digital Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Mammography Image Storage	1.2.840.10008.5.1.4.1.1.1.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Magnetic Resonance Image Storage	1.2.840.10008.5.1.4.1.1.4	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Nuclear Medicine Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.5	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Positron Emission	1.2.840.10008.5.1.4.1.1.128	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Tomography Image Storage		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Ultrasound Image Storage	1.2.840.10008.5.1.4.1.1.6.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Ultrasound Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.6	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Ultrasound Multi-frame Image Storage	1.2.840.10008.5.1.4.1.1.3.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2.4.91	SCP	None
Ultrasound Multi-frame Image Storage (Retired)	1.2.840.10008.5.1.4.1.1.3	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2.1	SCP	None
		JPEG baseline	1.2.840.10008.1.2.4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2.4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2.4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2.4.90	SCP	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2. 4.91	SCP	None
X-Ray Angiographic Bi- Plane Image Storage	1.2.840.10008.5. 1.4.1.1.12.3	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2. 1	SCP	None
		JPEG baseline	1.2.840.10008.1.2. 4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2. 4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2. 4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2. 4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2. 4.91	SCP	None
X-Ray Angiographic Image Storage	1.2.840.10008.5. 1.4.1.1.12.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2. 1	SCP	None
		JPEG baseline	1.2.840.10008.1.2. 4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2. 4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2. 4.70	SCP	None
		JPEG 2000 Lossless	1.2.840.10008.1.2. 4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2. 4.91	SCP	None
X-Ray Radiofluoroscopic Image Storage	1.2.840.10008.5. 1.4.1.1.12.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2. 1	SCP	None
		JPEG baseline	1.2.840.10008.1.2. 4.50	SCP	None
		JPEG extended	1.2.840.10008.1.2. 4.51	SCP	None
		JPEG lossless	1.2.840.10008.1.2. 4.70	SCP	None

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
		JPEG 2000 Lossless	1.2.840.10008.1.2. 4.90	SCP	None
		JPEG 2000 Lossless/lossy	1.2.840.10008.1.2. 4.91	SCP	None
Greyscale Softcopy Presentation State Storage	1.2.840.10008.5. 1.4.1.1.11.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2. 1	SCP	None
Key Image Note	1.2.840.10008.5. 1.4.1.1.88.59	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR, Little Endian	1.2.840.10008.1.2. 1	SCP	None
Mammography CAD SR IOD	1.2.840.10008.5. 1.4.1.1.88.50	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None

3.1.3.1.2.1 SOP Specific Conformance for SOP Class Verification

The *dcserver* Application Entity conforms to the DICOM Verification Service Class as an SCP.

3.1.3.1.2.2 SOP Specific Conformance for SOP Class Storage

The *dcserver* Application Entity conforms to the DICOM Storage Service Class as an SCP for the Abstract Syntax's listed in the table in section 3.1.3.1.2 at conformance level 2.

When a C-STORE operation is successful, the data has been stored in the database. The data is accessed through the workstation. The operator of the workstation determines the storage duration of the data.

The received images are grouped and sorted into Patient/Study/Series folders depending on values in the image. This is not directly related to the DICOM behaviour and is therefore not specified in this document.

In the case where the database is full, a status of 0xC001 is returned to the Storage SCU and the Information Object is discarded. The recovery action is to provide more storage space.

In the case where the image already exists in the workstation database (same Image instance UID), a status of 0xD000 is returned to the Storage SCU and the Information Object is discarded. A different returned value can be specified in a configuration file. The behaviour can apply to all associations or to specific Application entities.

If the Information Object Instance does not match any accepted Abstract Syntax, a status code of 0xA800 is returned. Recovery consists of altering the configuration of the remote or local Application Entity.

The attribute (0000,0902) contains a descriptive message to explain error returns.

The *dcserver* performs some validation before storing the image. Failure of a validation results in the return of status 0xC001 in the C-STORE response message. The following validations are performed:

- Invalid or missing orientation vector values (0020,0037) results in rejection of MR and CT Information Objects.
- A missing Photo Interpretation attribute (0028,0004) results in rejection of any modality Information Object.
- Missing Bits Allocated, Bits Used and High Bit attributes result in rejection of any modality Information Object.
- Any attributes with a value longer than what DICOM specifies are rejected.

The storage implementation performs the following coercions:

- If Pixel Padding Value (0028,0120) is present, the pixel values are adjusted accordingly.

The workstation supports a subset of the features that the Greyscale Softcopy Presentation State IOD supports. Namely, it supports the Presentation State, Display Shutter, Overlay Plane, Displayed Area, Graphic Annotation, Spatial Transformation, Graphic Layer, Modality LUT, and Softcopy VOI LUT modules.

3.1.3.1.3 Presentation Context Acceptance Criterion

The *dcserver* accepts Storage SOP Class Presentation Contexts if they are configured in the Application Entity configuration file. The possible Presentation Contexts are listed in section 3.1.3.1.2.

3.1.3.1.4 Transfer Syntax Selection Policies

The *dcserver* presently supports the default DICOM Little-endian Transfer Syntax.

3.1.3.2 Query Association Request

3.1.3.2.1 Associated Real-World Activity

The *dcserver* searches the attached database for the requested Information Objects described in the C_FIND identifier and returns a response for each match.

3.1.3.2.2 Presentation Context Table

The following table lists the possible Presentation Contexts. The Application Entity configuration file specifies which of these Presentation Contexts are actually used in a specific configuration.

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Patient Root Query / Retrieve Model C-FIND	1.2.840.10008.5.1.4.1.2.1.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
Study Root Query / Retrieve Model C-FIND	1.2.840.10008.5.1.4.1.2.2.1	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None

3.1.3.2.2.1 SOP Specific Conformance for Patient Root Query/Retrieve Model - FIND

The *dcserver* Application Entity conforms to the DICOM Patient Root Query/Retrieve Service Class as an SCP for the Abstract Syntax's listed in the table in section 3.1.3.2.2. The table in section 3.1.2.1.2.1 defines the accepted search keys.

A response is returned for each match found in the attached database.

Possible response status values are:

Refused	Out of resources	A700
Failed	Identifier does not match SOP Class	A900
	Unable to Process	C000
Cancel	Terminated due to Cancel Request	FE00
Success	Matching completed	0000
Pending	Matches are continuing	FF00

If the Information Object Instance does not match any accepted Abstract Syntax, a status code of 0xA800 is returned. Recovery consists of altering the configuration of the remote or local Application Entity.

The attribute (0000,0902) contains a descriptive message to explain error returns.

3.1.3.2.2.2 SOP Specific Conformance for Study Root Query/Retrieve Model - FIND

The *dcserver* Application Entity conforms to the DICOM Study Root Query/Retrieve Service Class as an SCP for the Abstract Syntax's listed in the table in section 3.1.3.2.2. The table in section 3.1.2.1.2.2 defines the accepted search keys.

A response is returned for each match found in the attached database.

Possible response status values are:

Refused	Out of resources	A700
Failed	Identifier does not match SOP Class	A900

	Unable to Process	C000
Cancel	Terminated due to Cancel Request	FE00
Success	Matching completed	0000
Pending	Matches are continuing	FF00

If the Information Object Instance does not match any accepted Abstract Syntax, a status code of 0xA800 is returned. Recovery consists of altering the configuration of the remote or local Application Entity.

The attribute (0000,0902) contains a descriptive message to explain error returns.

3.1.3.2.3 Presentation Context Acceptance Criterion

The *dcserver* accepts SOP Class contexts if they are configured in the Application Entity configuration file. The possible Presentation Contexts are listed in section 3.1.3.1.2.

3.1.3.2.4 Transfer Syntax Selection Policies

The *dcserver* presently supports the default DICOM Implicit Little-endian Transfer Syntax.

3.1.3.3 Move Association Request

3.1.3.3.1 Associated Real-World Activity

The *dcserver* initiates an association to the destination Application Entity specified in the C_MOVE command message. The *dcserver* then extracts the requested Information Objects described in the C_MOVE identifier from the attached database and performs C_STORE operations on the destination association.

3.1.3.3.2 Presentation Context Table

The following table lists the possible Presentation Contexts. The Application Entity configuration file specifies which of these Presentation Contexts are actually used in a specific configuration.

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Patient Root Query / Retrieve Model - MOVE	1.2.840.10008.5.1.4.1.2.1.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None
Study Root Query / Retrieve Model - MOVE	1.2.840.10008.5.1.4.1.2.2.2	Implicit VR, Little Endian	1.2.840.10008.1.2	SCP	None

3.1.3.3.2.1 SOP Specific Conformance for Patient Root Query/Retrieve Model - MOVE

The *dcserver* Application Entity conforms to the DICOM Patient Root Query/Retrieve Service Class as an SCP for the Abstract Syntax's listed in the table in section 3.1.3.3.2.

A response is returned for each Information Object sent to the destination Application Entity.

Possible response status values are:

Refused	Out of resources	A700
	Unable to perform sub-operations	A702
	Move Destination Unknown	A801
Failed	Identifier does not match SOP Class	A900
	Unable to Process	C000
Cancel	Terminated due to Cancel Request	FE00
Success	sub-operations completed	0000
Warning	sub-operations completed, 1 or more failures	B000
Pending	Matches are continuing	FF00

The attribute (0000,0902) contains a descriptive message to explain error returns.

3.1.3.3.2 SOP Specific Conformance for Study Root Query/Retrieve Model - MOVE

The *dcserver* Application Entity conforms to the DICOM Study Root Query/Retrieve Service Class as an SCP for the Abstract Syntax's listed in the table in section 3.1.3.3.2.

A response is returned for each Information Object sent to the destination Application Entity.

Possible response status values are:

Refused	Out of resources	A700
	Unable to perform sub-operations	A702
	Move Destination Unknown	A801
Failed	Identifier does not match SOP Class	A900
	Unable to Process	C000
Cancel	Terminated due to Cancel Request	FE00
Success	sub-operations completed	0000
Warning	sub-operations completed, 1 or more failures	B000
Pending	Matches are continuing	FF00

The attribute (0000,0902) contains a descriptive message to explain error returns.

3.1.3.3.3 Presentation Context Acceptance Criterion

The *dcserver* accepts SOP Class contexts if they are configured in the Application Entity configuration file. The possible Presentation Contexts are listed in section 3.1.3.3.2.

3.1.3.3.4 Transfer Syntax Selection Policies

The *dcserver* presently supports the default DICOM Implicit Little-endian Transfer Syntax.

3.2 Hcserver Application Entity Specification

The *hcserver* represents a single Application Entity. It acts independently of other DICOM applications that may be running on the same system. The *hcserver* can support printing to multiple DICOM printers at the same time, each printer being uniquely identified by an Application Entity Title.

The *hcserver* provides standard conformance to the following DICOM 3.1 SOP Classes as an SCU:

SOP Class Name	SOP Class UID
Basic Greyscale Print Management Meta	1.2.840.10008.5.1.1.9

3.2.1 Association Establishment Policies

3.2.1.1 General

The *hcserver* maintains a separate association with each DICOM SCP. It releases the association with the DICOM SCP if no operation is done on the association in a selected time period.

3.2.1.2 Number of Associations

There is no limit on the number of associations maintained simultaneously with one or different DICOM SCPs.

3.2.1.3 Asynchronous Nature

The application does not support asynchronous operations and will not perform asynchronous window negotiation.

3.2.1.4 Implementation Identifying Information

The *hcserver* implementation class UID is 2.16.124.113531.1.1.1.

The *hcserver* implementation version name is ISG_HCS_V1.0.96.

3.2.2 Association Initiation Policy

The *hcserver* maintains a list of valid print devices and can present that list to the applications upon request. When the application submits a print job designated for a listed print device to the *hcserver*, the *hcserver* will request an association with the selected print device.

3.2.2.1 Print to remote printer

3.2.2.1.1 Associated Real World Activity

The application's print request causes the *hcserver* to initiate an association.

3.2.2.1.2 Proposed Presentation Contexts

The *hcserver* will propose one of the presentation contexts listed in the Presentation Context Table.

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name	UID		
Basic Greyscale Print Management	1.2.840.10008.5.1.1.9	Implicit VR, Little Endian	1.2.840.10008.1.2	SCU	None

3.2.2.1.2.1 SOP Specific Conformance to Basic Greyscale Print Management Meta SOP Class

The *hcs*server supports the following mandatory SOP classes that are defined under the Basic Greyscale Print Management Meta SOP Class:

SOP Class Name	SOP Class UID
Basic Film Session	1.2.840.10008.5.1.1.1
Basic Film Box	1.2.840.10008.5.1.1.2
Basic Greyscale Image Box	1.2.840.10008.5.1.1.4
Printer	1.2.840.10008.5.1.1.16

The *hcs*server supports the following optional SOP class attributes and DIMSE services for the Basic Greyscale Print Management Meta SOP Class.

SOP Class	DIMSE Service	Optional Attribute	Tag
Basic Film Session SOP Class	N-CREATE	Number of Copies	(2000,0010)
		Print Priority	(2000,0020)
		Medium Type	(2000,0030)
		Film Destination	(2000,0040)
		Film Session Label	(2000,0050)
		Memory Allocation	(2000,0060)
Basic Film Box SOP Class	N-CREATE	Film Orientation	(2010,0040)
		Film Size ID	(2010,0050)
		Magnification Type	(2010,0060)
		Max Density	(2010,0130)
		Configuration Information	(2010,0150)
		Smoothing Type	(2010,0080)
		Border Density	(2010,0100)
		Empty Image Density	(2010,0110)
Min Density	(2010,0120)		

SOP Class	DIMSE Service	Optional Attribute	Tag
		Trim	(2010,0140)
	N-DELETE		
Basic Greyscale Image Box SOP Class	N-SET	Polarity	(2020,0020)
		Magnification type	(2010,0060)
		Smoothing type	(2010,0080)
		Requested Image Size	(2020,0030)
Printer SOP Class	N-GET		

3.2.2.1.2.1 Basic Film Session SOP Class attributes

The hserver supports the following mandatory and optional attribute values in this SOP class:

Attribute Name	Tag	Supported values
Number of Copies	(2000,0010)	Integer string
Print Priority	(2000,0020)	HIGH,MED,LOW
Medium Type	(2000,0030)	PAPER,CLEAR FILM,BLUE FILM
Film Destination	(2000,0040)	MAGAZINE, PROCESSOR
Film Session Label	(2000,0050)	Long string
Memory Allocation	(2000,0060)	Integer string

3.2.2.1.2.1.2 Basic Film Box SOP Class attributes

The hserver supports the following mandatory and optional attribute values in this SOP class:

Attribute Name	Tag	Supported values
Image Display Format	(2010,0010)	STANDARD, ROW, COL, SLIDE, SUPERSLIDE, CUSTOM
Film Orientation	(2010,0040)	PORTRAIT, LANDSCAPE
Film Size ID	(2010,0050)	8INX10IN, 10INX14IN, 14INX14IN, 24CMX24CM, 10INX12IN, 11INX14IN, 14INX17IN, 24CMX30CM
Magnification Type	(2010,0060)	REPLICATE, BILINEAR, CUBIC, NONE
Smoothing Type	(2010,0080)	SCP specific
Border Density	(2010,0100)	BLACK, WHITE, i where i represents the desired density in hundredths of OD
Empty Image Density	(2010,0110)	BLACK, WHITE, i where i represents the desired density in hundredths of OD
Min Density	(2010,0120)	Unsigned short
Max Density	(2010,0130)	Unsigned short
Trim	(2010,0140)	YES, NO
Configuration Information	(2010,0150)	SCP specific

3.2.2.1.2.1.3 Basic Greyscale Image Box SOP Class attributes

The hcserver supports the following mandatory and optional attribute values in this SOP class:

Attribute Name	Tag	Supported values
Image Position	(2020,0010)	Unsigned short
Polarity	(2020,0020)	NORMAL, REVERSE
Magnification Type	(2010,0060)	REPLICATE, BILINEAR, CUBIC, NONE
Smoothing Type	(2010,0080)	SCP specific
Requested Image Size	(2020,0030)	Unsigned short
Preformatted Greyscale Image Sequence	(2020,0110)	
>Samples per Pixel	(0028,0002)	1
>Photometric Interpretation	(0028,0004)	MONOCHROME1, MONOCHROME2
>Planar configuration	(0028,0006)	1
>Rows	(0028,0010)	Unsigned short
>Columns	(0028,0011)	Unsigned short
>Pixel Aspect Ratio	(0028,0034)	1:1
>Bits Allocated	(0028,0100)	8
>Bits Stored	(0028,0101)	8
>High Bit	(0028,0102)	7
>Pixel Representation	(0028,0103)	000H(unsigned integer)
>Pixel Data	(7FE0,0010)	Other Byte String

3.2.2.1.2.1.4 Printer SOP Class attributes

The hcserver makes use of the following attributes and attributes values in this SOP class:

Attribute Name	Tag	Supported values
Printer Status	(2110,0010)	NORMAL, WARNING FAILURE
Printer Status Info	(2110,0020)	SUPPLY EMPTY, SUPPLY LOW, RECEIVER FULL, FILM JAM
Printer Name	(2110,0030)	Long string
Manufacturer	(0008,0070)	Long string
Manufacturer Model Name	(0008,1090)	Long string
Device Serial Number	(0018,1000)	Long string
Software Versions	(0018,1020)	Long string(s)

3.2.3 Association Acceptance Policy

The *hcserver* does not accept associations.

4. Communication Profiles

4.1 Supported Communication Stacks (Parts 8,9)

The workstation DICOM services provide DICOM V3.0 TCP/IP Network Communication Support as defined in Part 8 of the DICOM standard.

4.1.1 TCP/IP Stack

The workstation DICOM services inherit its TCP/IP stack from the Windows OS system upon which they execute.

4.1.1.1 API

The implementation uses Berkeley style sockets.

4.1.1.2 Physical Media Support

The implementation is not dependent on the physical medium used for the TCP/IP network.

5. Extensions/Specialization's/Privatization's

5.1 Standard/Extended/Specialized/Private SOPs

Not applicable

5.2 Private Transfer Syntax's

No Private Transfer Syntax's are used.

6. Configuration

6.1 Dcserver Application Entity Configuration

The Query/Retrieve and Storage SOP Classes to accept are configurable, globally or Application Entity Title specific.

The Query/Retrieve and Storage SOP Classes to propose are configurable, globally or Application Entity Title specific.

The Transfer Syntax's are configurable for each SOP Class, globally or SOP Class specific.

A configuration parameter is supplied to control matching of Calling Application Entity Title to a value in the configuration file.

A configuration parameter is supplied to control matching of Called Application Entity Title to a value in the configuration file.

A configuration parameter is supplied to allow Application Entity Title specific association related tracing output to be created for connection troubleshooting.

A configuration parameter is supplied to allow Application Entity Title specific DIMSE tracing output to be created for message troubleshooting.

Application entity host names can be specified as either IP address or host name.

The number of associations that can be initiated is configurable.

The number of associations that can be accepted is configurable.

The port number to listen on for association requests is configurable.

Mapping between attributes in DICOM Information Objects and the target database is runtime configurable.

6.2 Hcserver Application Entity Configuration

Application entity host names can be specified as either IP address or host name.

The destination printer host name and port number is configurable. Multiple printers can be configured.

The film layout formats supported per printer are configurable.

The association timeout per printer is configurable.

7. Support of Extended Character Sets

This implementation supports the following extended character set:

ISO-IR 100 = Latin alphabet No. 1, supplementary set.